



Oxygen simplified

EverFlo oxygen concentrator

The smaller, lighter EverFlo oxygen concentrator from Respironics is making a big impression on providers and patients. Its sleek, ergonomic design makes EverFlo easier to store, move, and carry than conventional concentrators. And it is one of the quietest and most energy efficient devices available today. With no filter for patients to maintain, a durable metal oxygen outlet, and fewer moving parts than other concentrators, EverFlo is designed to deliver consistent, trouble-free operation.

PHILIPS

RESPIRONICS



Molded handle that's easy to pick up and transport with one or both hands



Humidifier bottle holder compatible with all bottle styles




Internal tubing storage keeps tubing conveniently available for use



Metal cannula connector for greater durability



Recessed flow meter reduces the risk of breakage

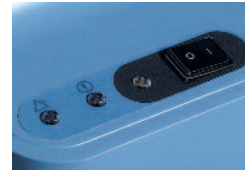

Interface panel with on/off switch and alarm indicators

Product specifications			
	Standard version	Low power version	IKK version
Power	230V/50Hz	230V/50Hz	230V/50Hz
Liter flow	0.5-5 liters per min	0.5-5 liters per min	0.5-5 liters per min
Sound level	43 dBA	43 dBA	< 40 dBA
Oxygen purity	90% - 96%	87% - 96%	87% - 96%
Weight	14 kg (32 lbs)	14 kg (32 lbs)	15 kg (33 lbs)
Dimensions (HxWxD)	580mm (22.8"), 380mm (15"), 240mm (9.5")	580mm (22.8"), 380mm (15"), 240mm (9.5")	580mm (22.8"), 380mm (15"), 240mm (9.5")
OPI (oxygen percentage indicator) alarm levels	Low oxygen: 82%, very low oxygen: < 70%	Low oxygen: 82%, very low oxygen: < 70%	Low oxygen: 82%, very low oxygen: < 70%
Required maintenance	Compressor intake filter: 2 yrs. Micro disk filter: as needed	Compressor intake filter: 2 yrs. Micro disk filter: as needed	Compressor intake filter: 2 yrs. Micro disk filter: as needed
Power consumption	< 300 W	Up to 3 lpm: 285 W; < 300 W at 5 lpm	< 300 W
Part numbers	1020006 Europe 1020009 South Africa 1020011 Italy 1020010 Australia 1020012 Argentina 1020013 South Korea 1020016 Hong Kong 1020017 Switzerland	1020008 UK, Ireland	1020007 Germany

Respironics and EverFlo are trademarks of Respironics, Inc and its affiliates. All rights reserved.



©2009 Koninklijke Philips Electronics N.V. All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

CAUTION: US federal law restricts these devices to sale by or on the order of a physician.

Broudy KB 00/00/09 MCI 4102636 PN 1039011

Philips Healthcare is part of Royal Philips Electronics

Respironics Asia Pacific
+65 6298 1088
Respironics Australia
+61 (2) 9666 4444
Respironics Europe, Middle East, Africa
+33 1 47 52 30 00
Respironics United Kingdom
+44 800 1300 845

Philips Respironics
1010 Murry Ridge Lane
Murrysville, PA 15668

Customer Service
+1 724 387 4000
+1 800 345 6443 (toll free, US only)
www.philips.com/respironics